Lifting armchairs

| 1. Introduction | page 504 |
|---|----------|
| 2. Protection elements and supports | page 506 |
| 3. Choosing the upholstery color | page 507 |
| 4. Dual motor models | page 508 |
| 5. Four motor models | page 514 |
| 6. Extra large model | page 516 |
| 7. Manual relax model | page 517 |
| 8. The Back School advices | page 519 |

RealeCollection Introduction

"When we have a moment to relax ourselves, we should avoid sinking into a too soft couch, but rather use an anatomic armchair which may allow us to abandon ourselves and relax our mind. Our spine will benefit greatly since disc pressure is reduced and the physiological shape is maintained." Prof. Benedetto Toso - Toso Back School Program (Italy)

Sturdiness is guaranteed by the use of high quality mechanisms and motoring, properly tested and standardized. With poplar or beech frames and plywood, while the armrests are made of particle boards. All screws have thread-locking compound to prevent unscrewing of the parts while moving.

The padding is made of high density foam rubber which prevents the change of shape over time. The thickness of the seat and the backrest are different enabling a measured support for the whole body.

The cover is made of artificial leather or microfiber. All covers are heavy-wear resistant and the upholstery finishings are high-quality.

The wheels enable easy transferring of the armchair on any surface and are highly resistant to static and dynamic stress.

Personalization is possible with several accessories. Furthermore, some changes to the measure of the structure, to the thickness of the padding, and the upholstery accessories can be applied providing they are feasible.

Safety is guaranteed by the use of motors with a 28V power adapter in accordance with the EU Directive 89/366. The adapter is placed outside the armchair and this reduces the exposure to electromagnetic fields. Last, some armchairs are certified as medical devices too.

Maximun safety load may vary according to the different models from a minimum of 110 kg to a maximum of 300 kg, with the Extra Large version.

The four-motor version is available for some models to allow not only full independence of the movement of the backrest and the legrest but also to adjust the whole armchair in height and to reach basculamento and Trendelenburg positions. The maximum load is increased up to 150 kg.

ETING ARMCHAIRS









The mechanical and moving parts of the relax armchairs and aid to lifting of the Chinesport program undergo load and performing tests to ensure reliability and long lasting performance. The guarantee period is three years for all models.

RealeCollection Introduction

Made in Italy

ERGONOMIC CHARACTERISTICS

Height of the seat:

it allows the elderly to use less energy while getting up, as well as to keep a straight angle for the lower limbs.

Depth of the seat:

it avoids the person from sliding forward.

Height of the backrest:

it provides a good support at the upper and lower back regions as well as at cervical area while resting.

Height of the armrests:


it decreases the preassure on the lower-back spine discs and to losen the shoulders;

Tilt while lifting:

it allows to adopt an upright position with the least effor preventing sliding due to eccessive tilt or too much flexion and spinal cord level;

Position of the legrest in relation to the hip:

When the legrest is completely raised, the legs are higher than the hip and this creates a lesser pressure on the discs and an easier blood flow.

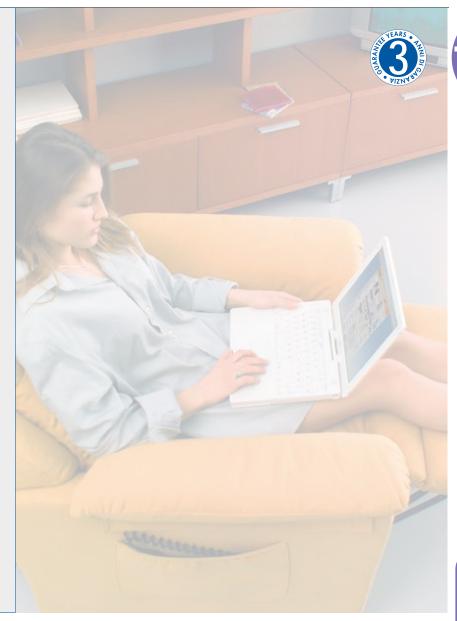

The backrest of some armchairs can be removed so they can be transferred easily..





DID YOU KNOW THAT...? Chinesport relax riser armchairs, with two or more motors make it easier to get up and sitting down for those who are impaiered in movement. Considering this important feature, as an aid, patients can get advantage of a reduced VAT percentage all along the Italian territory, when a Local Health Care System specialist issues a prescription, and when presenting an invalid certificate issued by an Official Government Institution.

RealeCollection Accessories

PROTECTION ELEMENTS AND SUPPORTS



01222. HEAD PROTECTION It is a microfiber cloth which protects the backrest. It is removable so as to be easily cleaned. It is always included in Lucy and July microfiber models.



01223. ARMREST PROTECTION They are microfiber covers to protect the armrest, easy to remove and to clean. Always included in Lucy and July microfiber models. Sold in pairs.



AC0752. **HAND CONTROL HOLDER** The hand control can be left at a fixed position during the lifting, and in special way at the completed raising stage. By this way the end user can help himself to get a standing position during the lifting with both arms on the armrests as lever.



01225. CERVICAL CUSHION It is a cushion for cervical support which can be adjusted in height. It can be easily removed and cleaned.



01226. CERVICAL U-SHAPE CUSHION This cervical support is an important restraint for the head on its sides when needed. It is padded and height adjustable.



01227. LOW-BACK CUSHION It is a cushion for the low-back support which can be chosen in the same color of the chair, as well as the other elements considering the upholstery.



POSTURAL CUSHIONS

Cushions of different sizes. They are elastic and can change shape, as they always return to the original shape, they are the ideal support for the low back (and lateral) area and the upper and lower limbs. The upholstery is light and washable.



01449 POSTURAL CUSHION 100 - ø cm 19 x 17 01450 POSTURAL CUSHION 110 - ø cm 19 x 31 01451 POSTURAL CUSHION 120 - ø cm 25 x 37 01452 POSTURAL CUSHION 8 - cm 40 x 20 x 13 h 01453 POSTURAL CUSHION 88 - cm 50 x 28 x 13 h

01426 POSTURAL SET The Postural kit consists of 5 cushions, one of each kind, as shown above, and a carry bag of cm 33 x 33 x 56 h size.

UPHOLSTERY CHOICE

Gi tto Colection **ARTIFICIAL LEATHER Standard Line** Giotto is the flame retardant fabric made for residential COMFORT and contract upholstery. Its resistant structure and AND ELASTICITY MAINTENANCE the easy cleaning make Giotto particularly suitable for hospitals, nursing homes, spas, rehabilitation centers HEAVY-WEAR UV-RADIATION and gyms. Giotto has a leather-effect surface with a RESISTANCE RESISTANCE brilliant finishing, in a wide range of colors all available on stock. Giotto contains Vinyzene, an anti-microbial WATERPROOF **RESISTANCE TO** agent added in the first steps of the production MOLD AND BACTERIA process. This creates a defensive barrier against mites, bacteria and fungus, avoiding smells and allergies. **W7 W8** WK WS WA WN DARK BROWN WHITE **IVORY** GREY WB **W4** WT W1 W6 WE DARK GREEN LIGHT GREEN CREAM WZ WG WF WH W9 WO ORANGE CORAL RED BURGUNDY **FUCHSIA** LILAC WR **W**3 WL WM WP **W2** LIGHT CORNFLOWER BLUE LIGHT BLUE MARINE BLUE DARK BLUE ANTHRACITE BLACK

F MICROFIBER

Tekna Plus handwash or in a washing machine up to 40°C. Do not wring, do not centrifugate. Let it dry inside out, do not expose to direct sunlight. Iron, inside out at low temperatures. For spots (oil, fat, cosmetics, ink, coffee, liquors, chewing-gum) it is recommended to use ethyl alcohol), (colorless) tamper with a detergent neutral solvent and wash carefully.



02047.? LUCY

Dual motor riser and relax armchair. It has the possibility to adjust the backrest and legrest independently, by means of a low-voltage (28V) hand control. This chair helps the user to adopt a reading posture, for watching television, relaxing or resting. It also rises and leans forward for aiding the user while sitting down or getting up. The person can adopt any intermediate position while getting up. The chair can be chosen with artificial leather or cloth upholstery. Whenever the choice is cloth, covers for the armrests and backrest are included. The covers are removable and easy to clean. Three-year guarantee period.

02048.? LUCY R

This model has the same characteristics as code 02047 but it also includes casters for transferring. Pedal control for activating the casters. Three-year guarantee period.



This model has a removable backrest, for easy transferring.







Detail of the pedal control for activating the casters of the armchair - Lucy R model code 02048.

| TECHNICAL SPECIFIC | CATION |
|---|-----------------------|
| CODE | |
| Backrest adjustment | electric, independent |
| Legrest adjustment | electric, independent |
| Riser function | electric |
| Transferring casters activation | Pedal control |
| Height of the seat from ground (cm) | 46 |
| Width of seat (cm) | 50 |
| Depth of seat (cm) | 50 |
| Height of the seat when raising a person (cm) | 56 |
| Armchair overall width (cm) | 75 |
| Armchair overall depth when closed (cm) | 84 |
| Armchair overall depth when open (cm) | 160 |
| Armchair height (cm) | 105 |
| Maximum load (kg) | 110 |
| Weight of armchair (kg) | 57 |

ACCESSORIES

Protection elements:

01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION AC0752.? CONTROL HOLDER

Low-back/cervical supports: 01225.? CERVICAL CUSHION 01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION

Postural cushion:

| ? COLOR CHOICE Please, specify the code for the chosen upholstery and color, as well as the armchair code on your purchase order. | |
|--|--|
| ARTIFICIAL LEATHER | |

| WA | WN | W8 | W7 | WK | ws | | W4 | WT | W1 | W6 | WE | WZ | WG | WF | WH | W9 | WQ | WR | W2 | W3 | WL | wм | WP |
|-------|-------|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| MICRO | FIBER | | | | | | | | | | | | | | | | | | | | | | |
| F8 | FC | FS | FD | F4 | FZ | FG | FH | FF | FQ | FT | F6 | F1 | FR | F2 | F3 | FL | FP | | | | | | |

CE

110 Kg

02090.? JULY

Dual motor riser and relax chair. It can be adjusted by means of a low-voltage hand control. This chair helps the user to adopt a reading posture, for watching television, relaxing or resting. It also rises and leans forward for aiding the user while sitting down or getting up. Whenever the choice is cloth, covers for the armrests and backrest are included. Armrests can be removed to allow the user to comfortably move sideways from/to a wheelchair. Three-year guarantee period.

02089.? JULY R


This model has the same characteristics as code 02090 but it also includes casters for transferring. Pedal control for activating the casters. Three-year guarantee period.



This model has a removable backrest, for easy transferring.

Detail of removable armrest.



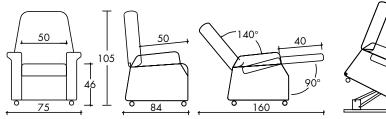



The July armchair features removable armrests on both sides for greater accessibility.

| IECHNICAL SPECIF | ICATION |
|---|-------------------------|
| CODE | |
| Backrest adjustment | electric, independent |
| Legrest adjustment | elettrica, independente |
| Riser function | electric |
| Transferring casters activation | Pedal control |
| Height of the seat from ground (cm) | 46 |
| Width of seat (cm) | 50 |
| Depth of seat (cm) | 50 |
| Height of the seat when raising a person (cm) | 56 |
| Armchair overall width (cm) | 75 |
| Armchair overall depth when closed (cm) | 84 |
| Armchair overall depth when open (cm) | 160 |
| Armchair height (cm) | 105 |
| Maximum load (kg) | 110 |
| Weight of armchair (kg) | 57 |

ACCESSORIES

Protection elements: 01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION AC0752.? CONTROL HOLDER

Low-back/cervical supports: 01225.? CERVICAL CUSHION 01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION

Postural cushion: 01449 POSTURAL CUSHION 100 01450 POSTURAL CUSHION 110 01451 POSTURAL CUSHION 120 01452 POSTURAL CUSHION 8 01453 POSTURAL CUSHION 88 01426 POSTURAL SET


? COLOR CHOICE Please, specify the code for the chosen upholstery and color, as well as the armchair code on your purchase order.

| ARTIFI | CIAL LE | ather | | | | | | | | | | | | | | | | | | | | | |
|------------|---------|-------|--------------|------|-----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| WA | WN | | W7 | wк | ws | | W4 | WT | W1 | W6 | WE | wz | WG | WF | WH | W9 | WQ | WR | W2 | W3 | WL | WM | WP |
| MICRO | FIBER | | | | | | | | | | | | | | | | | | | | | | |
| F 8 | FC | FS | FD | F4 | FZ | FG | FH | FF | FQ | FT | F6 | F1 | FR | F2 | F3 | FL | FP | | | | | | |
| | | | Harry | 1.62 | 1.5 | | | | | | | | | | | | | | | | | | |



01230.? FANNY 2EL

Dual motor riser and relax chair. It has the possibility to adjust the backrest and legrest independently, by means of a low-voltage (28V) hand control. This chair helps the user to adopt a reading posture, for watching television, and for better digesting. It also rises and leans forward for aiding the user while sitting down or getting up. The person can adopt any intermediate position while getting up. Lift is carried out safely thanks to the sturdy frame that rests on the floor and due to a smooth and continuous movement as well. This version allows the backrest to be even more reclined, and the leg rest to be lengthened by further 10 cm in relax and digesting position, it is possible to tilt the chair to raise the legs and improve the blood flow and decrease the pressure on discs. It is provided with a large pockets on the right side and with four swivel casters. The seat and leg rest are made from the same cushion and cover, which ensures a steady and uniform support for the lower limbs. Three-year guarantee period.

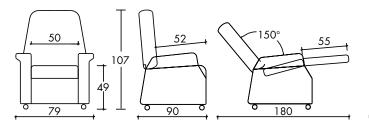



Π

120 Kg



| TECHNICAL SPECIFICA | ATION |
|---|-----------------------|
| CODE | |
| Backrest adjustment | electric, independent |
| Legrest adjustment | electric, independent |
| Riser function | electric |
| Transferring casters activation | Pedal control |
| Height of the seat from ground (cm) | 49 |
| Width of seat (cm) | 50 |
| Depth of seat (cm) | 52 |
| Height of the seat when raising a person (cm) | 68 |
| Armchair overall width (cm) | 79 |
| Armchair overall depth when closed (cm) | 90 |
| Armchair overall depth when open (cm) | 180 |
| Armchair height (cm) | 107 |
| Maximum load (kg) | 110 |
| Weight of armchair (kg) | 52 |

ACCESSORIES

Protection elements: 01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION AC0752.? CONTROL HOLDER

Low-back/cervical supports: 01225.? CERVICAL CUSHION 01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION

Postural cushion:

| 2 COLOR CHOICE Please, specify the code for the chosen upholstery and color, as well as the armchair code on your pu | urchase order. |
|--|----------------|
| | |

| AKII | | LEAI | HEK | | | | | | | | | | | | | | | | | | | | | |
|------|-------|------|-----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| WA | V | WN | W8 | W7 | WK | ws | | W4 | WT | W1 | W6 | WE | WZ | WG | WF | WH | W9 | WQ | WR | W2 | W3 | WL | wм | WP |
| MICR | OFIBE | ER | | | | | | | | | | | | | | | | | | | | | | |
| F8 | | FC | FS | FD | F4 | FZ | FG | FH | FF | FQ | FT | F6 | F1 | FR | F2 | F3 | FL | FP | | | | | | |

Made in Italy

01234.? ANDREA 2EL

Dual motor riser and relax chair. It has the possibility to adjust the backrest and legrest independently, by means of a low-voltage (28V) hand control. This chair helps the user to adopt a reading posture, for watching television, and for better digesting. It also rises and leans forward for aiding the user while sitting down or getting up. The person can adopt any intermediate position while getting up. Lift is carried out safely thanks to the sturdy frame that rests on the floor and due to a smooth and continuous movement as well. This version allows the backrest to be even more reclined, and the leg rest to be lengthened by further 10 cm in relax and digesting position, it is possible to tilt the chair to raise the legs and improve the blood flow and decrease the pressure on discs. It is provided with two side pockets and with four swivel casters. The seat and leg rest are made from the same cushion and cover, which ensures a steady and uniform support for the lower limbs. All the cushions (backrest, armrests, seat and legrest) can be removed, their cover can also be removed and cleaned in accordance with the technical characteristics of upholstery materials. Three-year guarantee period.





| | • <u>• 185</u> • | |
|---|-----------------------|-----------|
| TECHNICAL SPECIFIC | ATION | ACCESS |
| CODE | | Protectio |
| Backrest adjustment | electric, independent | 01222.? |
| Legrest adjustment | electric, independent | 01223.? |
| Riser function | electric | AC0752.3 |
| Transferring casters activation | Pedal control | Low-bac |
| Height of the seat from ground (cm) | 50 | 01225.? (|
| Width of seat (cm) | 50 | 01226.? |
| Depth of seat (cm) | 50 | 01227.? |
| Height of the seat when raising a person (cm) | 68 | Postural |

| Width of seat (cm) | 50 |
|---|-----|
| Depth of seat (cm) | 50 |
| Height of the seat when raising a person (cm) | 68 |
| Armchair overall width (cm) | 90 |
| Armchair overall depth when closed (cm) | 100 |
| Armchair overall depth when open (cm) | 185 |
| Armchair height (cm) | 112 |
| Maximum load (kg) | 110 |

65

110 Kg



68

SORIES on elements:

HEAD PROTECTION ARMREST PROTECTION **?** CONTROL HOLDER

k/cervical supports: **CERVICAL CUSHION** U-SHAPE CUSHION LOW-BACK CUSHION cushion:

01449 POSTURAL CUSHION 100 01450 POSTURAL CUSHION 110 01451 POSTURAL CUSHION 120 01452 POSTURAL CUSHION 8 01453 POSTURAL CUSHION 88 01426 POSTURAL SET

COLOR CHOICE Please, specify the code for the chosen upholstery and color, as well as the armchair code on your purchase order.

Weight of armchair (kg)

50

112

50


| ARTIFIC | CIAL LEAT | THER | | | | | | | | | | | | | | | | | | | | | |
|---------|-----------|------|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| WA | WN | | W7 | WK | ws | | W4 | WT | W1 | W6 | WE | WZ | WG | WF | WH | W9 | WQ | WR | W2 | W3 | WL | wм | WP |
| MICRO | FIBER | | | | | | | | | | | | | | | | | | | | | | |
| F8 | FC | FS | FD | F4 | FZ | FG | FH | FF | FQ | FT | F6 | F1 | FR | F2 | F3 | FL | FP | | | | | | |

01832.? PETRA 2EL

Dual motor riser and relax chair. It has the possibility to adjust the backrest and legrest independently, by means of a low-voltage (28V) hand control. This chair helps the user to adopt a reading posture, for watching television, and for relaxing. It also rises and leans forward for aiding the user while sitting down or getting up. The person can adopt any intermediate position while getting up. All the backrest cushions, armrests, seat and legrest covers can be removed and cleaned in accordance with the technical characteristics of upholstery materials. Three-year guarantee period.





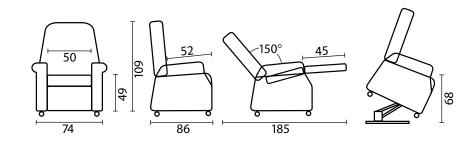



| TECHNICAL SPECIFICATION | | | | | | | |
|---|-----------------------|--|--|--|--|--|--|
| CODE | | | | | | | |
| Backrest adjustment | electric, independent | | | | | | |
| Legrest adjustment | electric, independent | | | | | | |
| Riser function | electric | | | | | | |
| Transferring casters activation | Pedal control | | | | | | |
| Height of the seat from ground (cm) | 49 | | | | | | |
| Width of seat (cm) | 50 | | | | | | |
| Depth of seat (cm) | 52 | | | | | | |
| Height of the seat when raising a person (cm) | 68 | | | | | | |
| Armchair overall width (cm) | 75 | | | | | | |
| Armchair overall depth when closed (cm) | 86 | | | | | | |
| Armchair overall depth when open (cm) | 185 | | | | | | |
| Armchair height (cm) | 109 | | | | | | |
| Maximum load (kg) | 110 | | | | | | |
| Weight of armchair (kg) | 58 | | | | | | |

ACCESSORIES

Protection elements: 01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION AC0752.? CONTROL HOLDER

Low-back/cervical supports: 01225.? CERVICAL CUSHION 01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION


Postural cushion:

| COLOR CHOICE Please, specify the code for the chosen upholstery and color, as well as the armchair code on your purchase order. | |
|--|--|
| | |

| A | KIIFIC | IAL LEA | INEK | | | | | | | | | | | | | | | | | | | | | |
|---|--------|---------|------|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| | WA | WN | W8 | W7 | WK | ws | | W4 | WT | W1 | W6 | WE | WZ | WG | WF | WH | W9 | WQ | WR | W2 | W3 | WL | WM | WP |
| N | ICRO | FIBER | | | | | | | | | | | | | | | | | | | | | | |
| | F8 | FC | FS | FD | F4 | FZ | FG | FH | FF | FQ | FT | F6 | F1 | FR | F2 | F3 | FL | FP | | | | | | |





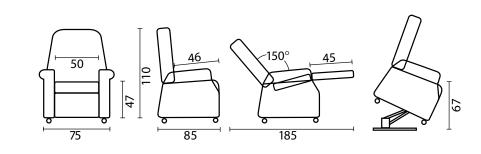



| TECHNICAL SPECIF | ICATION |
|---|-----------------------|
| CODE | 01834 |
| Backrest adjustment | electric, independent |
| Legrest adjustment | electric, independent |
| Riser function | electric |
| Transferring casters activation | Pedal control |
| Height of the seat from ground (cm) | 47 |
| Width of seat (cm) | 50 |
| Depth of seat (cm) | 46 |
| Height of the seat when raising a person (cm) | 67 |
| Armchair overall width (cm) | 75 |
| Armchair overall depth when closed (cm) | 85 |
| Armchair overall depth when open (cm) | 185 |
| Armchair height (cm) | 110 |
| Maximum load (kg) | 110 |
| Weight of armchair (kg) | 56 |

ACCESSORIES Protection elements: 01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION

AC0752.? CONTROL HOLDER Low-back/cervical supports: 01225.? CERVICAL CUSHION

01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION Postural cushion:

01449 POSTURAL CUSHION 100 01450 POSTURAL CUSHION 110 01451 POSTURAL CUSHION 120 01452 POSTURAL CUSHION 8 01453 POSTURAL CUSHION 88 01426 POSTURAL SET

? COLOR CHOICE Please, specify the code for the chosen upholstery and color, as well as the armchair code on your purchase order.

| ARTIFICIAL LE | ATHER | | | | | | | | | | | | | | | | | | | | | |
|---------------|-------|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| WA WN | W8 | W7 | wк | ws | | W4 | WT | W1 | W6 | WE | WZ | WG | WF | WH | W9 | WQ | WR | W2 | W3 | WL | wм | WP |
| MICROFIBER | | | | | | | | | | | | | | | | | | | | | | |
| F8 FC | FS | FD | F4 | FZ | FG | FH | FF | FQ | FT | F6 | F1 | FR | F2 | F3 | FL | FP | | | | | | |

01833. PETRA 4EL-BASCULA

Four-motor riser and relax chair. It provides complete independent movements for backrest, footrest and seat angle as to the floor. These movements are electrically adjusted by a low voltage (28V) hand control. This allows the chair to recline at the correct posture for reading and watching TV, or for getting up, while the four motors also make it possible to rise, tilt and adopt the Trendelenburg position. All the backrest cushions, armrests, seat and legrest covers can be removed and cleaned in accordance with the technical characteristics of upholstery materials. Three-year guarantee period.





| | | 150° 45 185 | |
|--|--|----------------|--|
|--|--|----------------|--|

| TECHNICAL SPECIFIC | ATION |
|---|-----------------------|
| CODE | |
| Backrest adjustment | electric, independent |
| Legrest adjustment | electric, independent |
| Riser function | electric |
| Transferring casters activation | Pedal control |
| Height of the seat from ground (cm) | 49 |
| Width of seat (cm) | 50 |
| Depth of seat (cm) | 52 |
| Height of the seat when raising a person (cm) | 68 |
| Armchair overall width (cm) | 74 |
| Armchair overall depth when closed (cm) | 86 |
| Armchair overall depth when open (cm) | 185 |
| Armchair height (cm) | 109 |
| Maximum load (kg) | 150 |
| Weight of armchair (kg) | 65 |

ACCESSORIES

Protection elements: 01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION AC0752.? CONTROL HOLDER

Low-back/cervical supports: 01225.? CERVICAL CUSHION 01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION

Postural cushion:

| COLOR CHOICE Please, specify the code for the c | | | | |
|--|---------------------|---------------------------|---------------------|----------------------|
| | haaan | مطدمه المتترمة بتمامه امت | annachair cada an i | wave average and a |
| E CUTUR CHUICE Please specify the code for the c | nosen hinnoisierv (| no color as well as the | ormenoir code on | vour nurchose order |
| | | ind color, as woll as ino | | your porchase oraon. |

| ARTIFIC | ial leat | HER | | | | | | | | | | | | | | | | | | | | | |
|---------|----------|-----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| WA | WN | | W7 | wк | ws | | | WT | W1 | W6 | WE | wz | WG | WF | ₩Н | W9 | WQ | WR | W2 | W3 | WL | wм | WP |
| MICRO | FIBER | | | | | | | | | | | | | | | | | | | | | | |
| F8 | FC | FS | FD | F4 | FZ | FG | FH | FF | FQ | FT | F6 | F1 | FR | F2 | F3 | FL | FP | | | | | | |




CE

T 150 Kg

01835.? SALLY 4EL-BASCULA

Four-motor riser and relax chair. It provides complete independent movements for backrest, footrest and seat angle as to the floor. These movements are electrically adjusted by a low voltage (28V) hand control. This allows the chair to recline at the correct posture for reading and watching TV, or for getting up, while the four motors also make it possible for whole armchair to rise, tilt and adopt the Trendelenburg position. All the backrest cushions, armrests, seat and legrest covers can be removed and cleaned in accordance with the technical characteristics of upholstery materials. Three-year guarantee period.

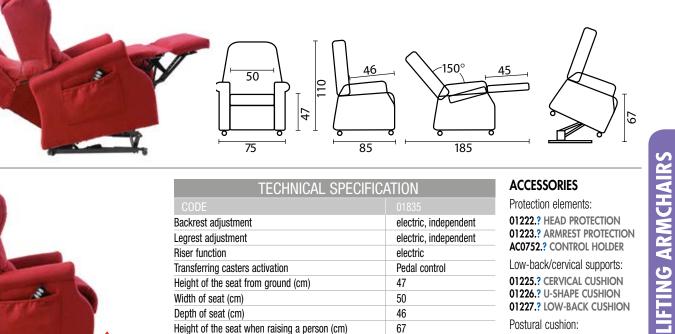



| TECHNICAL SPECIFIC | AHON |
|---|-----------------------|
| CODE | |
| Backrest adjustment | electric, independent |
| Legrest adjustment | electric, independent |
| Riser function | electric |
| Transferring casters activation | Pedal control |
| Height of the seat from ground (cm) | 47 |
| Width of seat (cm) | 50 |
| Depth of seat (cm) | 46 |
| Height of the seat when raising a person (cm) | 67 |
| Armchair overall width (cm) | 75 |
| Armchair overall depth when closed (cm) | 85 |
| Armchair overall depth when open (cm) | 185 |
| Armchair height (cm) | 110 |
| Maximum load (kg) | 150 |
| Weight of armchair (kg) | 63 |

ACCESSORIES Protection elements:

01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION AC0752.? CONTROL HOLDER

Low-back/cervical supports: 01225.? CERVICAL CUSHION 01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION Postural cushion:

01449 POSTURAL CUSHION 100 01450 POSTURAL CUSHION 110 01451 POSTURAL CUSHION 120 01452 POSTURAL CUSHION 8 01453 POSTURAL CUSHION 88 01426 POSTURAL SET

COLOR CHOICE Please, specify the code for the chosen upholstery and color, as well as the armchair code on your purchase order.

| ARTIFIC | ial leat | HER | | | | | | | | | | | | | | | | | | | | | |
|---------|----------|-----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| WA | WN | W8 | W7 | WK | ws | | W4 | WT | W1 | W6 | WE | wz | WG | WF | WH | W9 | WQ | WR | W2 | W3 | WL | wм | WP |
| MICRO | FIBER | | | | | | | | | | | | | | | | | | | | | | |
| F8 | FC | FS | FD | F4 | FZ | FG | FH | FF | FQ | FT | F6 | F1 | FR | F2 | F3 | FL | FP | | | | | | |

RealeCollection Four motor models - Extra large

01889. ? EXTRA LARGE

Four-motor riser and relax armchair of large size. It provides complete independent movements for backrest and footrest. These movements are electrically adjusted by a low voltage hand control. This allows the chair to recline at the correct posture for reading and watching TV, or relaxing, while the four motors also make it possible for whole armchair to rise, tilt and adopt the Trendelenburg position. Provided with a side pocket. The armchair rests on a steel supporting frame to ensure maximum strength to stand a load up to 300 kg. It is a medical device. Three-year guarantee period.





| (and a | |
|--------|---|
| | |
| | + |

| TECHNICAL SPECIF | ICATION |
|---|-----------------------|
| CODE | |
| Backrest adjustment | electric, independent |
| Legrest adjustment | electric, independent |
| Riser function | electric |
| Transferring casters activation | Pedal control |
| Height of the seat from ground (cm) | 52 |
| Width of seat (cm) | 65 |
| Depth of seat (cm) | 51 |
| Height of the seat when raising a person (cm) | 68 |
| Armchair overall width (cm) | 92 |
| Armchair overall depth when closed (cm) | 85 |
| Armchair overall depth when open (cm) | 185 |
| Armchair height (cm) | 114 |
| Maximum load (kg) | 300 |
| Weight of armchair (kg) | 66 |

114

52

<u>6</u>5

92

ACCESSORIES Protection elements:

160

48

01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION AC0752.? CONTROL HOLDER Low-back/cervical supports: 01225.? CERVICAL CUSHION 01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION

Γī 300 Kg

CE

68

Postural cushion:

| ? COLOR CHOICE Please, specify the code for the chosen upholstery and color, as well as the armo | hair code on your purchase order. |
|---|-----------------------------------|
|---|-----------------------------------|

| ART | FICI/ | al leat | HER | | | | | | | | | | | | | | | | | | | | | |
|-----|-------|---------|-----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| W | Ą | WN | | W7 | wк | ws | | | wт | W1 | W6 | WE | WZ | WG | WF | WH | W9 | WQ | WR | W2 | W3 | WL | WM | WP |
| MIC | ROFI | BER | | | | | | | | | | | | | | | | | | | | | | |
| F | 3 | FC | FS | FD | F4 | FZ | FG | FH | FF | FQ | FT | F6 | F1 | FR | F2 | F3 | FL | FP | | | | | | |

ealeCollection Manual models for relax

Made in Italy

02028. ? TRUDY


A rest armchair that is manually adjustable by exerting a light pressure on the armrests. In this armchair, the legrest can be lifted first and the backrest can be reclined at a later time. The backrest is balanced by unloading the weight onto the mechanism, so that the user can easily recline the back while keeping the position reached. It is therefore possible to reach a position for reading or relaxing. The backrest is contoured to provide support for a physiological lumbar lordosis, dorsal kyphosis and cervical lordosis. Three-year quarantee period.

02029. ? TRUDY 2

This model has the same characteristics as code 02028 but it also includes braked casters for transferring.





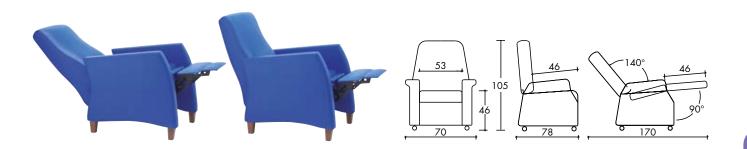



The backrest is removable which allows easy storing and transferring.

| TECHNICAL SPECIFICATION | | | | | |
|---|--------|--|--|--|--|
| CODE | | | | | |
| Backrest adjustment | manual | | | | |
| Legrest adjustment | manual | | | | |
| Riser function | nd | | | | |
| Transferring casters activation | nd | | | | |
| Height of the seat from ground (cm) | 46 | | | | |
| Width of seat (cm) | 53 | | | | |
| Depth of seat (cm) | 46 | | | | |
| Height of the seat when raising a person (cm) | nd | | | | |
| Armchair overall width (cm) | 70 | | | | |
| Armchair overall depth when closed (cm) | 78 | | | | |
| Armchair overall depth when open (cm) | 170 | | | | |
| Armchair height (cm) | 105 | | | | |
| Maximum load (kg) | 110 | | | | |
| Weight of armchair (kg) | 53 | | | | |

? COLOR CHOICE Please, specify the code for the chosen upholstery and color, as well as the armchair code on your purchase order.

| ARTIFICIAL LEATHER | | | | | | | | | | | | | | | | |
|--------------------|-------|-------|----|------|-------|-------|------|----|----|----|----|----|----|----|----|----|
| WA WN W8 W7 | WK WS | WB W4 | wт | W1 V | W6 WE | wz w | G WF | WH | W9 | WQ | WR | W2 | W3 | WL | WM | WP |
| MICROFIBER | | | | | | | | | | | | | | | | |
| F8 FC FS FD | F4 FZ | FG FH | FF | FQ | FT F6 | F1 FF | F2 | F3 | FL | FP | | | | | | |

ACCESSORIES Protection elements: 01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION AC0752.? CONTROL HOLDER

Postural cushion:

Low-back/cervical supports: 01225.? CERVICAL CUSHION 01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION

01449 POSTURAL CUSHION 100 01450 POSTURAL CUSHION 110

01451 POSTURAL CUSHION 120

01452 POSTURAL CUSHION 8 01453 POSTURAL CUSHION 88 01426 POSTURAL SET

RealeCollection Manual models for relax

33620.? FANNY M

A rest armchair that is manually adjustable by exerting a light pressure on the armrests. This chair helps the user to adopt a reading posture, for watching television, and for better digesting. Its width size has been considered so that it can be easily moved around the house. May the client make any special request, it can be customized. It has a right side pocket and four swivel non-marking rubber casters with brakes, for easy transfers. The seat and legrest are made from the same cushion and cover, which ensures a steady and even support for the lower limbs. It can be easily transported thanks to a slide-in system that allows for removing the backrest. Three-year guarantee period.




This model has a removable backrest, for easy transferring.





| TECHNICAL SPECIFI | ICATION |
|---|---------------|
| CODE | |
| Backrest adjustment | manual |
| Legrest adjustment | manual |
| Riser function | nd |
| Transferring casters activation | Pedal control |
| Height of the seat from ground (cm) | 47 |
| Width of seat (cm) | 50 |
| Depth of seat (cm) | 46 |
| Height of the seat when raising a person (cm) | nd |
| Armchair overall width (cm) | 75 |
| Armchair overall depth when closed (cm) | 85 |
| Armchair overall depth when open (cm) | 185 |
| Armchair height (cm) | 110 |
| Maximum load (kg) | 110 |
| Weight of armchair (kg) | 36 |



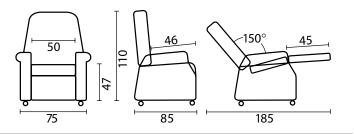

ACCESSORIES

Protection elements: 01222.? HEAD PROTECTION 01223.? ARMREST PROTECTION AC0752.? CONTROL HOLDER

Low-back/cervical supports: 01225.? CERVICAL CUSHION 01226.? U-SHAPE CUSHION 01227.? LOW-BACK CUSHION

Postural cushion:

| ? COLOR CHOICE Please, specify the code for the chosen upholstery | and color, as well as the armchair code on your purchase order. |
|--|---|
|--|---|

| ARTIFICIAL LEATHER | | | |
|--------------------|------------------|------------------|----------------------|
| WA WN W8 W7 WK WS | WB W4 WT W1 W6 W | E WZ WG WF WH W9 | WQ WR W2 W3 WL WM WP |
| MICROFIBER | | | |
| F8 FC FS FD F4 FZ | FG FH FF FQ FT F | 6 F1 FR F2 F3 FL | FP |

RealeCollection Advices from Back School


The 80% of vertebral ache is not caused by any specific pathology of the spinal cord, but by wrong movements and postures, psychological stress, unfit physical shape, overweight, obesity and smoking habits.

The Back School advices:

• always choose positions and movements in which there is a minor pressure on discs;

- keep physiological curves;
- change position frequently;



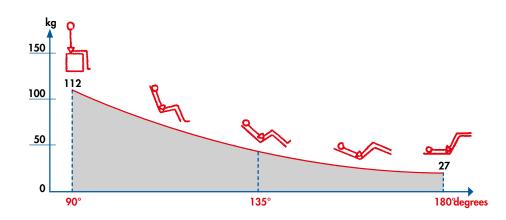

Alf L. Nachemson has demonstrated that the more reclined the backrest is, the lesser the pressure on the discs. To obtain the best results it is very important the backrest can keep the physiological curves.

The pictures show some examples of different activities which stress the low-back spinal cord. These are the causes of most low-back ache. It is better, then, to avoid repeating the same position when relaxing on the armchair.





